

# Pocket balances KERN TEE · TCB · CM · CM-C · TAB



## KERN TEE

Designer pocket balance at a bargain price



## KERN TCB

This is where the fun starts with pocket format weighing



## 1 KERN CM

Small, practical and protected in a hard case cover



## KERN TAB

High-precision pocket carat balance for jewellers and travelling salesmen

### STANDARD



### STANDARD



### OPTION



### STANDARD



### STANDARD



### OPTION



### STANDARD



### STANDARD



### OPTION



### STANDARD



### STANDARD



### OPTION



- Simple and convenient 2-key operation
- Cover made of shock proof plastic as protection against pressure and dust. Can also be used as a tare cup
- Note: Can be printed with your logo, from 100 pcs. Please ask for details

- Simple and convenient 2-key operation
- Small and handy
- Transport box and weighing pan included
- Revolving weighing plate (as overload protection)
- Note: Individual logo printing on request

- Hard case cover as protection against pressure and dust

### 1 CM:

- Integrated calculator

### 2 CM-C:

- Hard case cover as protection against pressure and dust
- Adjusting weight standard
- Draft shield against air movements standard
- Weighing pan standard


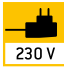

- Simple and convenient 4-key operation
- Can be switched over from g to ct, gn at the touch of a key
- Hard case cover as protection against pressure and dust
- Adjusting weight and weighing pan standard



2

KERN	TEE 150-1	TCB 200-1	CM 60-2N	CM 150-1N	CM 320-1N	CM 1K1N	CM 50-C2N	TAB 20-3	
<b>Weighing range [Max] g</b>	150 g	200 g	60 g	150 g	320 g	1000 g	10 g   50 ct	20 g   100 ct	
<b>Readout [d] g</b>	0,1 g	0,1 g	0,01 g	0,1 g	0,1 g	1 g	0,002 g   0,01 ct	0,001 g   0,005 ct	
<b>LCD display</b>	backlit, digit height 9 mm	Digit height 9 mm	Digit height 12 mm			Digit height 9 mm	backlit, digit height 9 mm		
<b>Dimensions of weighing plate</b>	WxD 60x64 mm	Ø 80 mm	WxD 70x80 mm			WxD 50x40 mm	WxD 36x42 mm		
<b>Overall dimensions</b>	WxDxH 67x100x22 mm	ØxH 80 mm	WxDxH 85x130x25 mm			WxDxH 95x133x33 mm			
<b>Power supply</b>	Batteries included, 2 x 1,5 V AAA	Batteries included, 2 x LR 44 AUTO-OFF function to preserve battery life, can be switched off	Batteries included, 2 x 1,5 V AAA AUTO-OFF function to preserve battery life, can be switched off			Batteries included, 4 x 1,5 V AAA operating time approx. 150 h			
			ca. 30 h	ca. 70 h	ca. 30 h				
<b>Net weight</b>	100 g	50 g	180 g			200 g			
<b>Permissible ambient temperature</b>	5 °C / 35 °C							10 °C / 30 °C	
<b>Option DAKKS Calibr. Certificate</b>	963-127								

# KERN Pictograms:

 <b>Internal adjusting:</b> Quick setting up of the balance's accuracy with internal adjusting weight (motordriven).	 <b>Piece counting:</b> Reference quantities selectable. Display can be switched from piece to weight.	 <b>Suspended weighing:</b> Load support with hook on the underside of the balance.
 <b>Adjusting program CAL:</b> For quick setting up of the balance's accuracy. External adjusting weight required.	 <b>Recipe level A:</b> Separate memory for the weight of the tare container and the recipe ingredients (net total).	 <b>Battery operation:</b> Ready for battery operation. The battery type is specified for each device.
 <b>Memory:</b> Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.	 <b>Recipe level B:</b> Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display.	 <b>Rechargeable battery pack:</b> Rechargeable set.
 <b>Alibi memory:</b> Electronic archiving of weighing results, complying with the 2009/23/EC standard.	 <b>Recipe level C:</b> Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display, adjustment of recipe when dosages are exceeded, multiplier function, barcode.	 <b>Universal mains adapter:</b> with universal input and optional input socket adapters for A) EU, GB B) EU, GB, CH, USA C) EU, GB, CH, USA, AUS
 <b>Data interface RS-232:</b> To connect the balance to a printer, PC or network.	 <b>Totalising level A:</b> The weights of similar items can be added together and the total can be printed out.	 <b>Mains adapter:</b> 230V/50Hz in standard version for EU. On request GB, USA or AUS version available.
 <b>RS-485 data interface:</b> To connect the balance to a printer, PC or other peripherals. High tolerance against electromagnetic disturbance.	 <b>Totalising level C:</b> Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display, adjustment of recipe when dosages are exceeded, multiplier function, barcode recognition.	 <b>Power supply:</b> Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, USA or AUS on request.
 <b>USB data interface:</b> To connect the balance to a printer, PC or other peripherals.	 <b>Weighing principle: Strain gauge</b> Electrical resistor on an elastic deforming body.	 <b>Weighing principle: Tuning fork</b> A resonating body is electromagnetically excited, causing it to oscillate.
 <b>Bluetooth* data interface:</b> To transfer data from the balance to a printer, PC or other peripherals.	 <b>Percentage determination:</b> Determining the deviation in % from the target value (100 %).	 <b>Weighing principle: Electromagnetic force compensation</b> Coil inside a permanent magnet. For the most accurate weighings.
 <b>WLAN data interface:</b> To transfer data from the balance to a printer, PC or other peripherals.	 <b>Weighing units:</b> Can be switched to e.g. non-metric units at the touch of a key. See balance model. Please refer to KERN's website for more details.	 <b>Weighing principle: Single cell technology</b> Advanced version of the force compensation principle with the highest level of precision.
 <b>Control outputs (optocoupler, digital I/O):</b> To connect relays, signal lamps, valves, etc.	 <b>Weighing with tolerance range:</b> Upper and lower limiting values can be programmed individually for e.g. dosing, sorting and portioning.	 <b>Verification possible:</b> The time required for verification is specified in the pictogram.
 <b>Interface for second balance:</b> For direct connection of a second balance.	 <b>Hold function:</b> (Animal weighing program) When the weighing conditions are unstable, a stable weight is calculated as an average value.	 <b>DAKkS calibration possible (DKD):</b> The time required for DAKkS calibration is shown in days in the pictogram.
 <b>Network interface:</b> For connecting the scale to an Ethernet network. With KERN products you can use a universal RS-232/LAN converter.	 <b>Protection against dust and water splashes IPxx:</b> The type of protection is shown in the pictogram.	 <b>Package shipment:</b> The time required for internal shipping preparations is shown in days in the pictogram.
 <b>Wireless data transfer:</b> between the weighing unit and the evaluation unit using an integrated radio module.	 <b>ATEX explosion protection:</b> Suitable for use in hazardous industrial environments, in which there is explosion danger. The ATEX marking is specified for each device.	 <b>Pallet shipment:</b> The time required for internal shipping preparations is shown in days in the pictogram.
 <b>GLP/ISO log:</b> The balance displays the weight, date and time, regardless of a printer connection.	 <b>Stainless steel:</b> The balance is protected against corrosion.	 <b>Warranty:</b> The warranty period is shown in the pictogram.
 <b>GLP/ISO log:</b> With weight, date and time. Only with KERN printers.		

## KERN – Precision is our business

To ensure the high precision of your balance KERN offers you the the appropriate test weight in the international OIML error limit classes E1-M3 from 1 mg - 2000 kg. In combination with a DAKkS calibration certificate the best pre-requisite for proper balance calibration.

The KERN DAKkS calibration laboratory today is one of the most modern and best-equipped DAKkS calibration laboratories for balances, test weights and force-measurement in Europe.

Thanks to the high level of automation, we can carry out DAKkS calibration of

balances, test weights and force-measuring devices 24 hours a day, 7 days a week.

### Range of services:

- DAKkS calibration of balances with a maximum load of up to 50 t
- DAKkS calibration of weights in the range of 1 mg – 2500 kg
- Database supported management of checking equipment and reminder service
- Calibration of force-measuring devices
- DAKkS calibration certificates in the following languages D, GB, F, I, E, NL, PL

## Your KERN specialist dealer: